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## CONTENTS

Features.....	15
Important safety instructions.....	16
Details description.....	18
Operating instruction.....	19
Working modes.....	20
Troubleshooting.....	22
Care & maintenance.....	23
Specifications.....	23

## FEATURES

- Deep kneading back massage
- 3 intensity levels of the vibrating massage for hips and buttocks
- Upper, lower or full back massage
- 4 rollers massage mechanism
- Timer: 15 minutes
- Removable cover helps to regulate the intensity of the massage
- For car, home and office

## IMPORTANT SAFETY INSTRUCTIONS

### **Please read this user's manual attentively.**

- Do not put hands or fingers near the massage mechanism while this product is running.
- Touching the mechanism can result in an injury from the rollers squeezing together.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use.
- If you feel any abnormal discomfort during the operation of this product, turn off the power immediately and do not use the product. Consult your physician.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Always keep the power cord away from heated surfaces.
- Comfortable clothing is recommended while using this product.
- Never drop or insert any object into any opening.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.
- Connect this appliance to a properly grounded outlet only.
- Do not use this product if noise above normal sound is heard.
- Do not stand on the appliance. Use only while seated.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- **Before the first use please remove the safety screw, located on the back side of the massage cushion with the use of a hexagonal key included.**
- Plug the massage cushion to the 12V connector only after you start the car.

**ATTENTION!**

- **Do not** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **Do not** operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- **Do not** sit on the controller.
- **Do not** pull by the AC cord.
- **Do not** carry this appliance by the power cord or use the cord as a handle.
- **Do not** reach for a product that has fallen into water. Unplug immediately.
- **Do not** use while bathing or in a shower.

**Persons who are undergoing medical treatment as well as the following persons should consult a physician before use:**


- Persons currently with lumbago, persons who have injured their legs, lower back, or neck in the past, persons who have numbness in their lower back, neck (persons with a condition, such as herniated disc, spondylolisthesis, herniated cervical vertebra, etc.).
- Persons having osteoarthritis, rheumatism, or gout.
- Persons having a bone condition, such as osteoporosis.
- Persons having a circulatory disorder (heart disease, angiopathy, hypertension, etc.).
- Persons having a respiratory organ disorder.

- Persons using a body implanted medical electronic device.
- Persons having malignant tumors.
- Persons having symptoms of a circulatory disorder (thrombosis, serious aneurysm, chronic varicosis, etc.), dermatitis, or infectious skin diseases.
- Persons having a sensory disorder due to an advanced peripheral circulatory disorder, such as that resulting from diabetes.
- Persons having skin injuries.
- Persons having a fever (38°C or higher) due to disease.
- Persons having an abnormal or crooked backbone.
- Persons who are pregnant or thought to be pregnant or who are menstruating.
- Persons who feel physical discomfort and require rest.
- Persons wishing to use the unit for rehabilitation purposes.
- Persons other than the above who do not feel physically normal.

Consult your physician  
before use!

# DETAILS DESCRIPTION
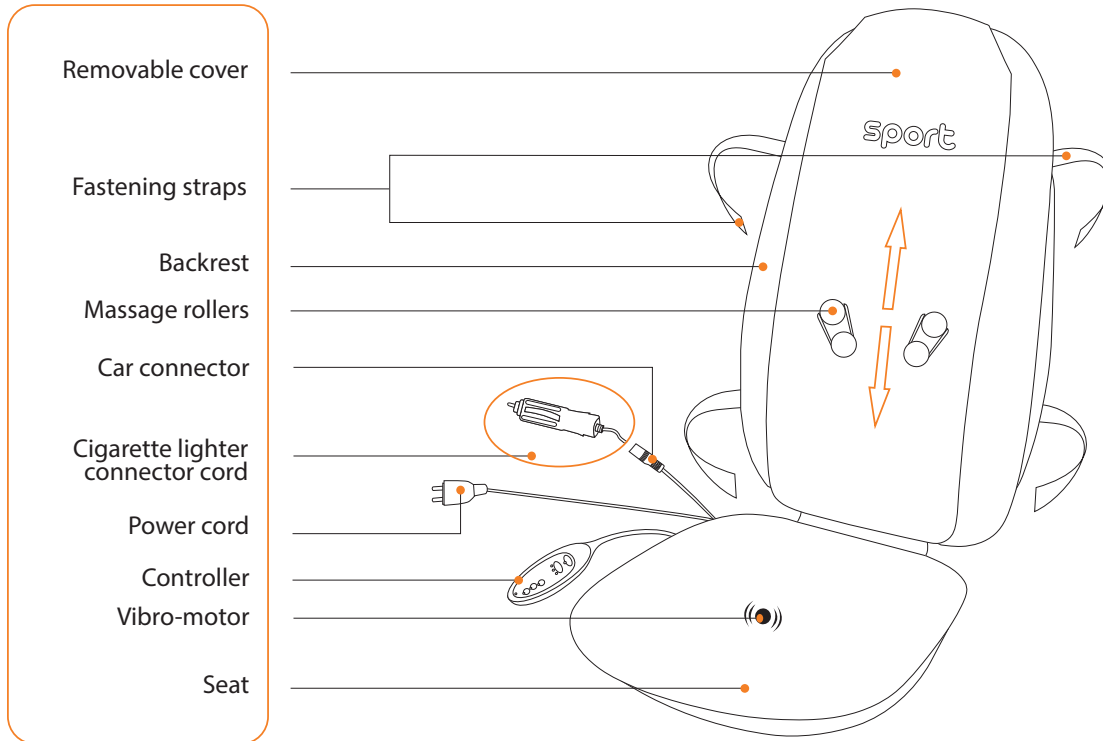
## OPERATING INSTRUCTION

### IN A CAR

### AT HOME

Remove the safety screw, located on the back side of the massage cushion with the use of a hexagonal key included. Slip the fastening straps over the seat or chair and adjust them to secure.


Plug one end of the cigarette lighter connector cord into the car connector on the massage cushion. Start the car and then plug the other end of the cigarette lighter connector cord into the 12V cigarette lighter connector in the car.

Plug the massage cushion into the 220V socket.

Press the **On/Off** button on the remote control.  
After finishing the massage, press the **On/Off** button on the remote control one more time.  
Unplug the power cord.

# WORKING MODES
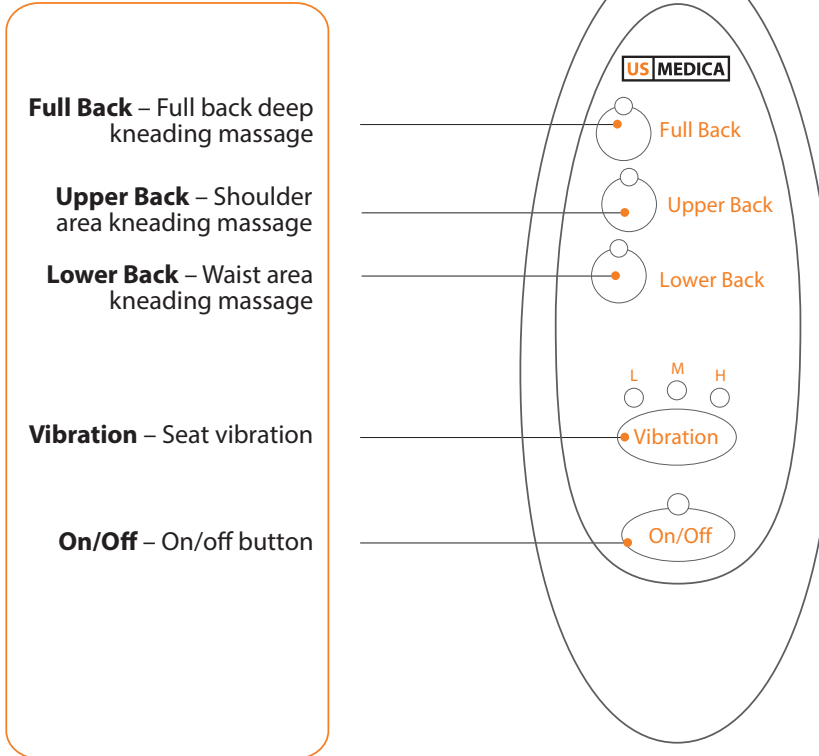
## Remote control

MODE	DESCRIPTION	BENEFITS
<p><b>BACK MASSAGE</b></p> 	<ul style="list-style-type: none"> <li>• Press <b>On/Off</b> button to turn on the massage cushion.</li> <li>• Press one of the three buttons to choose deep kneading back massage:</li> </ul> <p><b>Full Back</b> — Full back area massage  <b>Upper Back</b> — Upper back area massage  <b>Lower Back</b> — Lower back area massage</p>	<p>Relieves tension and soreness by lifting and stretching muscles, improves circulation and relaxes the back.</p>
<p><b>VIBRATION</b></p> 	<ul style="list-style-type: none"> <li>• Press <b>Vibration</b> button repeatedly to choose one of the three intensity levels of the seat vibration:</li> </ul> <p><b>L</b> — Low level of the vibration  <b>M</b> — Medium level of the vibration  <b>H</b> — High level of the vibration</p> <p><input checked="" type="radio"/> This mode can work individually and in a combination with back massage mode.</p>	<p>Improves circulation and helps to bring vital nutrients to the body. Helps to give temporary relief from minor muscle and joint pain and stiffness.</p>

## TROUBLESHOOTING

PROBLEM	CAUSE	SOLUTION
The device is very loud during massage	The massage cushion is not in full contact with the surface it was placed on	Ensure that the entire massage cushion is resting flat
Massage stopped	The preset maximum period of operation has elapsed	Let the device cool down for at least 15 min. before using again. Switch the device on again
	The connection to the mains has been interrupted	Check the plug connection and the contact of the main plug with the main supply
No massage	Appliance is not connected to the main power supply	Make all connections and switch on the device
	Device is not switched on	Turn the device on by pressing the <b>On/Off</b> button on the remote control

## CARE & MAINTENANCE

### Storage

- Do not wrap the power cord around the device before storing.
- Store the massage cushion in its original box and in a safe, dry and cool place.
- Store it out of the reach of children.
- Avoid contact with sharp edges or pointed objects as this might damage the upholstery.

### Cleaning

- Disconnect the device and allow it to cool down.
- Use a soft, damp cloth to clean the device.
- Use a soft brush to clean dirty connections and contacts.
- Dry the device with a soft cloth to minimise the potential water damage.
- Wait until the device is completely dry before using it again.

## SPECIFICATIONS

Name of the product	Massage cushion
Model	SPORT
Voltage	220 V, 12 V
Power	24 W

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US Medica LLC., 1201 N.Orange St, Wilmington, DE 19801, USA