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## CONTENTS

Features.....	17
Important safety instructions.....	18
Operating instruction .....	20
Working modes.....	22
Troubleshooting.....	25
Care & maintenance .....	26
Specifications.....	26

## FEATURES

- Deep kneading massage for back
- Rolling massage
- Vibrating massage for hips and buttocks
- Optional back heating
- Width function — adjust the rollers width to wide, mid or narrow
- Spot function — adjust the massage rollers' position in any particular spot of your back
- Use full back massage function, or just partial
- For car, home and office
- The removable cover helps to change the massage intensity
- Includes AC and DC adapters for use at home or in a car

## IMPORTANT SAFETY INSTRUCTIONS

### **Please read this user's manual attentively.**

- Do not put hands or fingers near the massage mechanism or rollers while this product is running.
- Touching the mechanism can result an injury from the rollers squeezing together.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
- The appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- If you feel any abnormal discomfort during the operation of this product, turn off the power immediately and do not use the product. Consult your physician.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Always keep the power cord away from heated surfaces.
- Never operate this product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- Comfortable clothing is recommended while using this product.
- Never drop or insert any object into any opening.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.
- This product was meant to be used in a relaxed position. Never force any body part into the backrest area while the rollers are moving.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.
- Connect this appliance to a properly grounded outlet only.
- Do not use this product if noise above normal sound is heard.
- Do not stand on or in appliance. Use only while seated.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Before use please remove the seal, located on the back side of the massage cushion with the use of a hexagonal key included.
- Plug the massage cushion to the 12V connector only after you start the car.

**SAVE THESE INSTRUCTIONS**

**ATTENTION!**

- **Do not** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **Do not** operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- **Do not** sit on the controller.
- **Do not** pull by the AC cord.
- **Do not** use outdoors.
- **Do not** carry this appliance by the power cord or use the cord as a handle.
- **Do not** reach for a product that has fallen into water. Unplug immediately.
- **Do not** use while bathing or in a shower.

**Persons who are undergoing medical treatment as well as the following persons should consult a physician before use:**


- Persons currently with lumbago, persons who have injured their legs, lower back, or neck in the past, persons who have numbness in their lower back, neck (persons with a condition, such as herniated disc, spondylolisthesis, herniated cervical vertebra, etc.).
- Persons having osteoarthritis, rheumatism, or gout.
- Persons having a bone condition, such as osteoporosis.
- Persons having a circulatory disorder (heart disease, angiopathy, hypertension, etc.).
- Persons having a respiratory organ disorder.
- Persons using a body implanted medical electronic device.
- Persons having malignant tumors.

- Persons having symptoms of a circulatory disorder (thrombosis, serious aneurysm, chronic varicosis, etc.), dermatitis, or infectious skin diseases.
- Persons having a sensory disorder due to an advanced peripheral circulatory disorder, such as that resulting from diabetes.
- Persons having skin injuries.
- Persons having a fever (38°C or higher) due to disease.
- Persons having an abnormal or crooked backbone.
- Persons who are pregnant or thought to be pregnant or who are menstruating.
- Persons who feel physical discomfort and require rest.
- Persons wishing to use the unit for rehabilitation purposes.
- Persons other than the above who do not feel physically normal.

Consult your physician  
before use!

# OPERATING INSTRUCTION
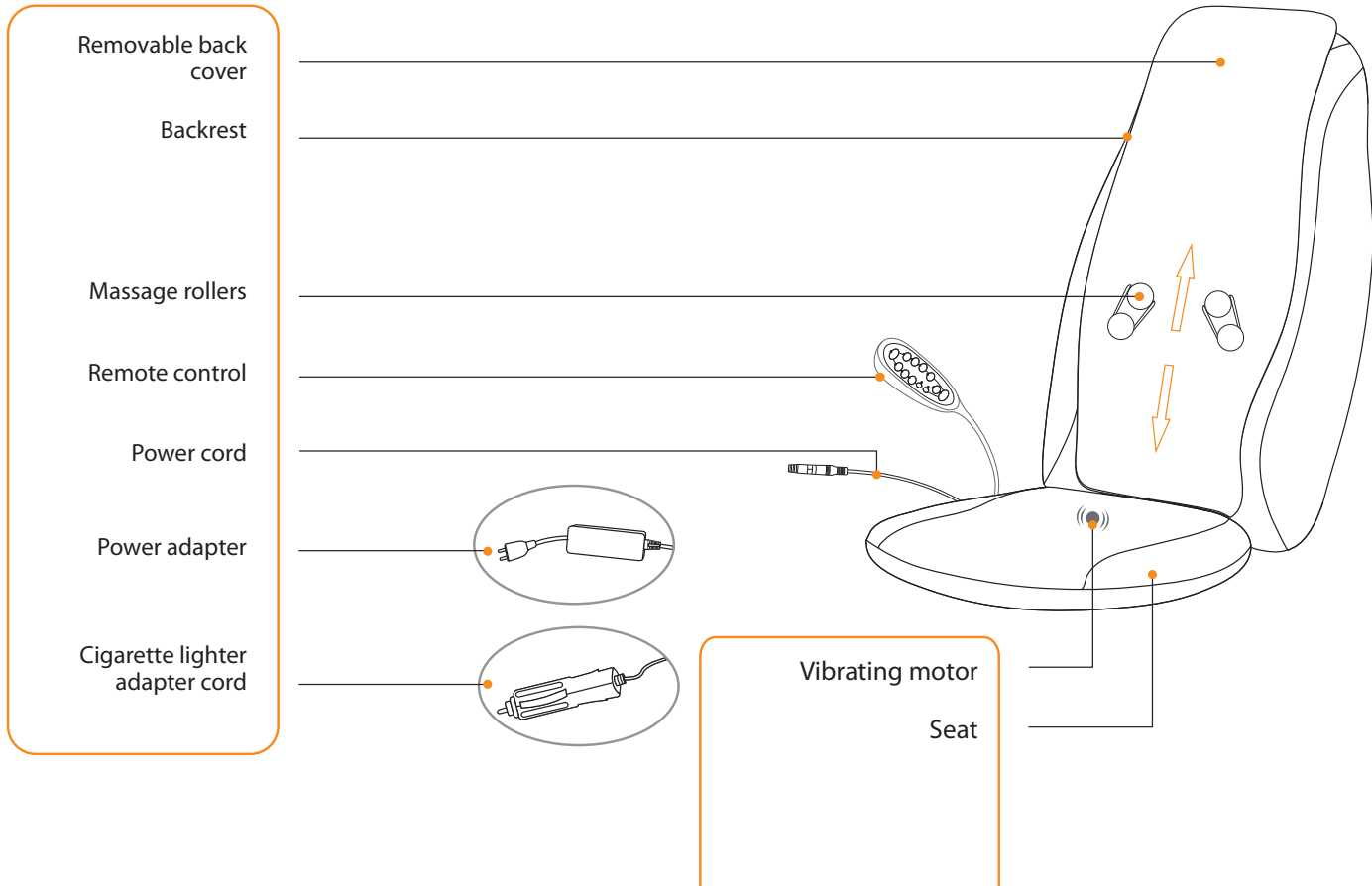
## IN A CAR

## AT HOME

Remove the seal, located on the back side of the massage cushion with the use of a hexagonal key included. Slip the fastening straps over the seat or chair and adjust them to secure.

Plug one end of the cigarette lighter connector cord into the power connector on the massage cushion. Start the car and then plug the other end of the cigarette lighter connector cord into the 12V cigarette lighter connector in the car.


Plug one end of the power adapter cord into the power cord connector on the massage cushion. Plug the other end into the 220V socket.

Press the power button on the remote control.  
After finishing the massage, press the power button on the remote control one more time.  
Unplug the power cord.

# WORKING MODES
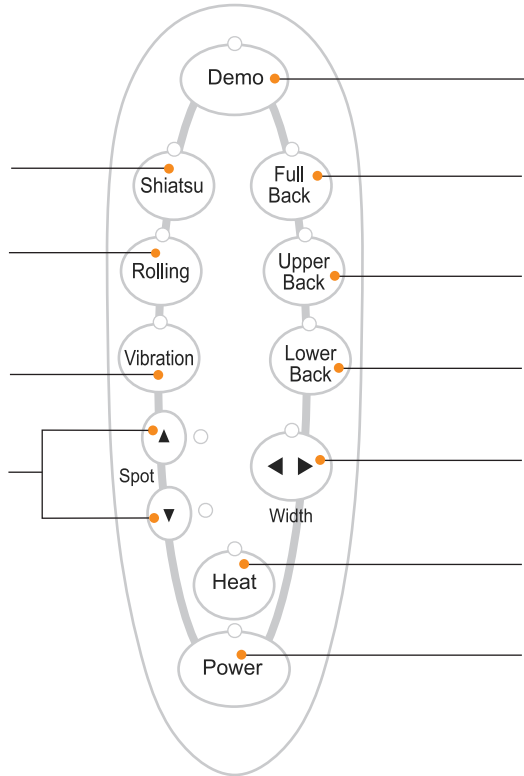
## Remote control

- Shiatsu** - Deep kneading mode
- Rolling** - Rolling mode
- Vibration** - Vibration mode
- Spot** - Particular spot massage

- Demo** - Demonstration mode
- Full Back** - Full back massage
- Upper Back** - Upper back area massage
- Lower Back** - Lower back area massage
- Width** - Adjust the position of the massage rollers
- Heat** - Heating function
- Power** - The power button

MODE	DESCRIPTION	BENEFITS
<p><b>DEMO</b></p> 	<p>Press <b>Demo</b> button for short demonstration of all the functions of the massage cushion.</p> <p>During the demonstration light indicators will be flashing on the remote control. After finishing the demonstration the massage cushion will automatically turn off.</p>	<p>Relieves tension and soreness by lifting and stretching muscles.</p>
<p><b>SHIATSU</b></p> 	<p>Press <b>Shiatsu</b> button for deep kneading your back on both sides of your spine simultaneously.</p> <p>During this mode you can adjust the position of the massage rollers in any particular spot of your back. Press and hold <b>Spot</b> «▲» or «▼» buttons to adjust the position:</p> <p>«▲» — the rollers move up;  «▼» — the rollers move down.</p> <p>⦿ Spot function is available only for <b>SHIATSU</b> mode.</p>	<p>Improves circulation and helps to bring vital nutrients to the spinal area.</p>
<p><b>ROLLING</b></p> 	<p>Press <b>Rolling</b> button and rollers will glide steadily along the curves of your back.</p> <p>During this mode you can adjust the rollers width to wide, mid or narrow. Press <b>Width</b> «↔» button to adjust the position.</p> <p>⦿ Width function is available only for <b>ROLLING</b> mode.</p>	<p>Warms and loosens your muscles, preparing your back for deeper massage.</p> <p>Relieves muscle tension, improves circulation and relaxes the back.</p>

MODE	DESCRIPTION	BENEFITS
<p><b>VIBRATION</b></p> 	<p>Press <b>Vibration</b> button for choosing one of three intensity levels of the vibration or for turning off the <b>VIBRATION</b> mode.</p> <p>☉ This mode can work individually and in a combination with <b>ROLLING</b> or <b>SHIATSU</b> mode.</p>	<p>Improves mobility and overall posture.</p>
<p><b>HEAT</b></p> 	<p>Heating function is available for <b>ROLLING</b> or <b>SHIATSU</b> mode. To turn on the function, during <b>ROLLING</b> or <b>SHIATSU</b>, press the <b>Heat</b> button. To turn off the function press the <b>Heat</b> button again.</p>	<p>Helps to give temporary relief from minor muscle and joint pain and stiffness.</p>

During **ROLLING** and **SHIATSU** modes you can choose the massage area by pressing one of the buttons:

**Full Back** – Full back massage

**Upper Back** – Upper back area massage

**Lower Back** – Lower back area massage

## TROUBLESHOOTING

PROBLEM	CAUSE	SOLUTION
The device is very loud during massage	The massage cushion is not fully In. The massage cushion is not in full contact with the surface it was placed on	Ensure that the entire massage cushion is resting flat
Massage stopped	The preset maximum period of operation has elapsed	Let the device cool down for at least 15 min before using again. Switch the device on again
	The connection to the mains has been interrupted	Check the plug connection and the contact of the main plug with the main supply
No massage	Appliance is not connected to the main power supply	Make all connections and switch on the device
	Device is not switched on	Turn the device on by pressing the <b>Power</b> button on the remote control

## CARE & MAINTENANCE

### Storage

- Do not wrap the power cord around the device before storing.
- Store the massage cushion in its original box and in a safe, dry and cool place.
- Store it out of the reach of children.
- Avoid contact with sharp edges or pointed objects as this might damage the fabric.

### Cleaning

- Disconnect the device and allow it to cool down.
- Use a soft, damp cloth to clean the device.
- Use a soft brush to clean dirty connections and contacts.
- Dry the device with a soft cloth to minimise potential water damage.
- Wait until the device is completely dry before using it again.

## SPECIFICATIONS

Name of the product	Massage cushion
Model	COMBO
Voltage	Entrance: 220 V, 12 V
Power	48 W
Includes AC and DC adapters for use at home or in a car	

[www.us-medica.com](http://www.us-medica.com)  
[www.us-medica.ru](http://www.us-medica.ru)

US Medica LLC., 1201 N.Orange St, Wilmington, DE 19801, USA